camry

CR 8155

(GB) user manual

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(LV) lietošanas instrukcija

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

(NL) handleiding

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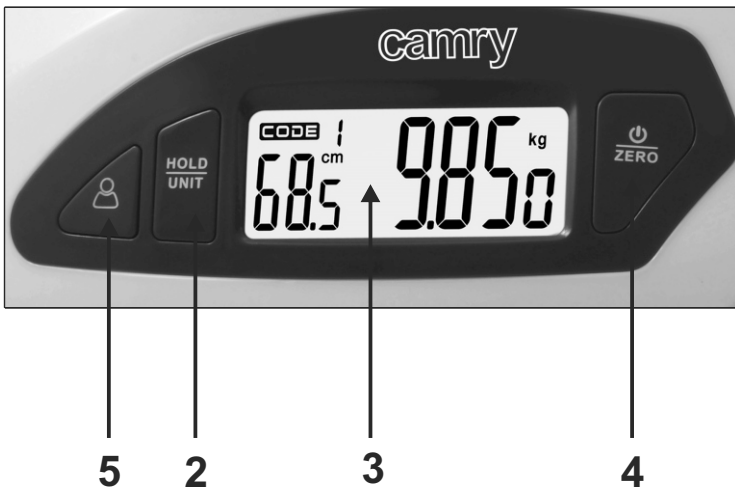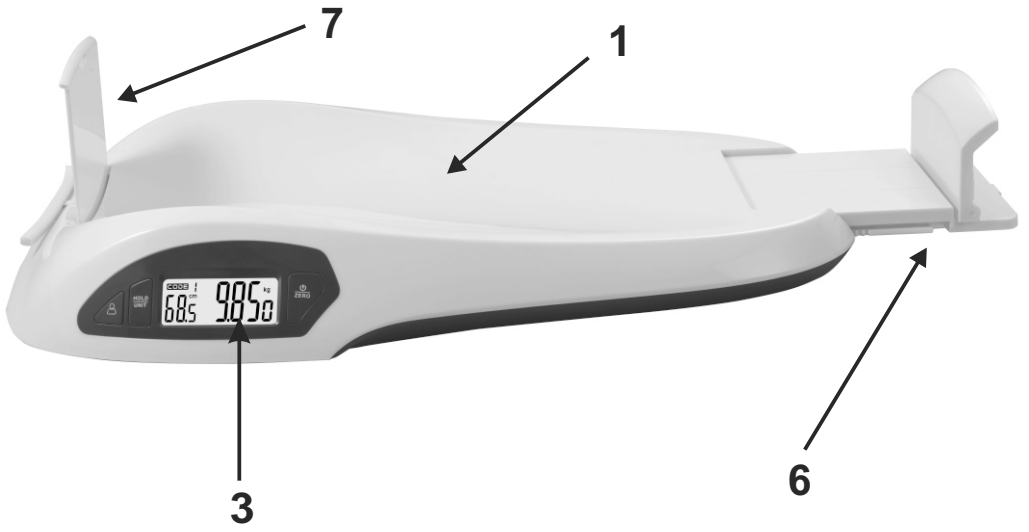SAFETY CONDITIONS. IMPORTANT INSTRUCTIONS ON SAFETY OF USE PLEASE READ CAREFULLY AND KEEP FOR FUTURE REFERENCE

The warranty conditions are different, if the device is used for commercial purpose.

1. Before using the product please read carefully and always comply with the following instructions. The manufacturer is not responsible for any damages due to any misuse.
2. The product is only to be used indoors. Do not use the product for any purpose that is not compatible with its application.
3. Please be cautious when using around children. Do not let the children to play with the product. Do not let children or people who do not know the device to use it without supervision.
4. **WARNING:** This device may be used by children over 8 years of age and persons with reduced physical, sensory or mental capabilities, or persons without experience or knowledge of the device, only under the supervision of a person responsible for their safety, or if they were instructed on the safe use of the device and are aware of the dangers associated with its operation. Children should not play with the device. Cleaning and maintenance of the device should not be carried out by children, unless they are over 8 years of age and these activities are carried out under supervision.
5. Never put the power cable, the plug or the whole device into the water. Never expose the product to the atmospheric conditions such as direct sun light or rain, etc.. Never use the product in humid conditions.
6. Never use the product with a damaged power cable or if it was dropped or damaged in any other way or if it does not work properly. Do not try to repair the defected product yourself because it can lead to electric shock. Always turn the damaged device to a professional service location in order to repair it. All the repairs can be done only by authorized service professionals. The repair that was done incorrectly can cause hazardous situations for the user.
7. Do not use the device near flammable materials.
8. Do not hold the device or the adapter when your hands are wet.
9. Batteries may leak if they are drained or were not used for an extended period of time. To protect the device and your health, replace them regularly and avoid skin contact with leaking batteries.
10. Use only AAA type manganese or alkaline batteries. Do not use other types of batteries. Do not use old and new batteries together. Immediately replace the old batteries with new ones.

DEVICE DESCRIPTION

- | | | |
|---------------------------|-------------------------------------|-------------------------|
| 1. Scale surface | 2. UNIT, HOLD button | 3. Display |
| 4. ON/OFF and ZERO button | 5. Button of the change of the user | 6. Frame of the measure |
| 7. Limiter of the measure | | |

REPLACING THE BATTERY

Remove the lid of the battery compartment at the bottom of the scale. Remove the old battery. Insert the new battery of the right type. Close the lid of the battery compartment.

BEFORE FIRST USE

Insert the correct batteries into the compartment (See "Replacing the battery")

SELECTING UNITS OF MEASURE

Press the ON/OFF button (4) to turn on the scale. Press the UNIT button (2) to change the units of measure. The unit will be visible on the display.

USING DEVICE

To ensure the accuracy of measurement, please use the scale on a fat and firm floor.

1. Press button "4" to turn on the scale, press button "5" to switch user and then long press again to confirm.

*Data will not be saved in guest mode (user no. "0").

2. Raise the bottom limiter to the purpose of the activation of the measurement of the height. LCD will show 46.0 the cm / 18.1 in as the initial value, if not, press the button "4" and next move the measuring plate in order to measure the height of the child.

3. Lay the baby precisely on the scale. The head should touch the upper limiter "7", however feet should touch the bottom limiter of

measurement board "6".

4. Press button "2" to hold the measurements of weight and height. Results flash twice for confirmation and are stored in scale memory. Comparison to previous data will be displayed in difference value repeatedly. In guest mode, measurements will be displayed for reference, but will not be saved. (Press "2" to cancel or let it turn off automatically without any operation.)

Memory of scale: Maximum 8 user supported. 15 records for each user. The first record will be overwritten by the latest record when the records are more than 15.

For weight measurement only, just skip step 2 & 3.

CLEANING AND MAINTENANCE

– Gently wipe the scale with a dry or slightly damp cloth.

MESSAGES

5 symbol – low battery

EEEE – excessive scale load

TECHNICAL DATA

Power: 4 x LR03 (AAA)

Capacity: 25 kg

Accuracy: 5g

Range Rulers: 46-80 cm

Accuracy: 1 mm

To protect your environment: please separate carton boxes and plastic bags and dispose them in corresponding waste bins. Used appliance should be delivered to the dedicated collecting points due to hazardous components, which may affect the environment. Do not dispose this appliance in the common waste bin.

Device is made in class II of insulation. Device is compliant with EU directives:

- Low voltage directive (LVD)

- Electromagnetic compatibility (EMC)

Device marked CE mark on rating label

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